



Booster Serum Product Introduction

REANZEN CO., LTD 2021.02.

CONTENTS


1. Intro
2. Key ingredients
3. Product Features
4. Product Line-up
5. How to use



About Us


Improve the Quality of Human life

- Bio-convergence & integration company
- Next-generation biotechnology research company focused on the core technology of convergence and integration with bio-materials, pharmaceuticals and devices
- Provide beauty solutions of innovative technologies with the knowhow of the founder and affiliated research laboratory
- Contribute to the happiness and healthy life of human beings, focusing on the medical esthetic business and derma care business

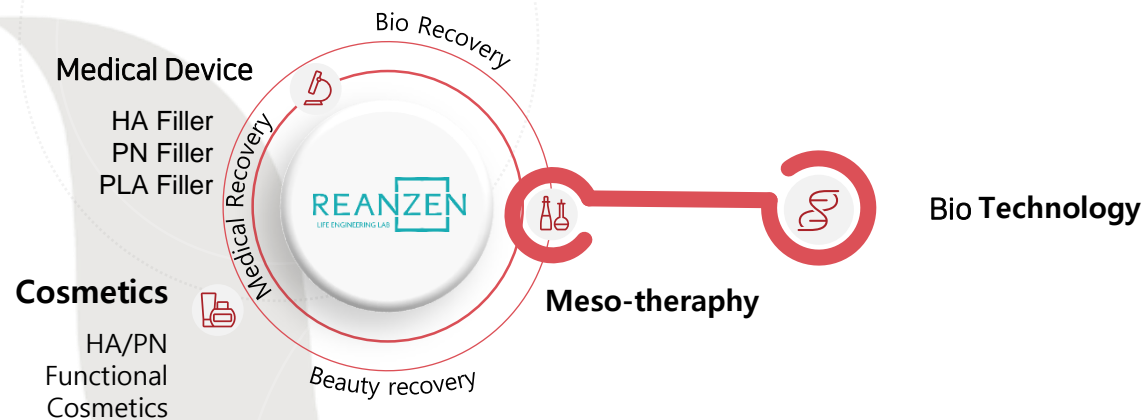

Bio-convergence & integration specialized company, REANZEN

REVS YOUTH LABORATORY

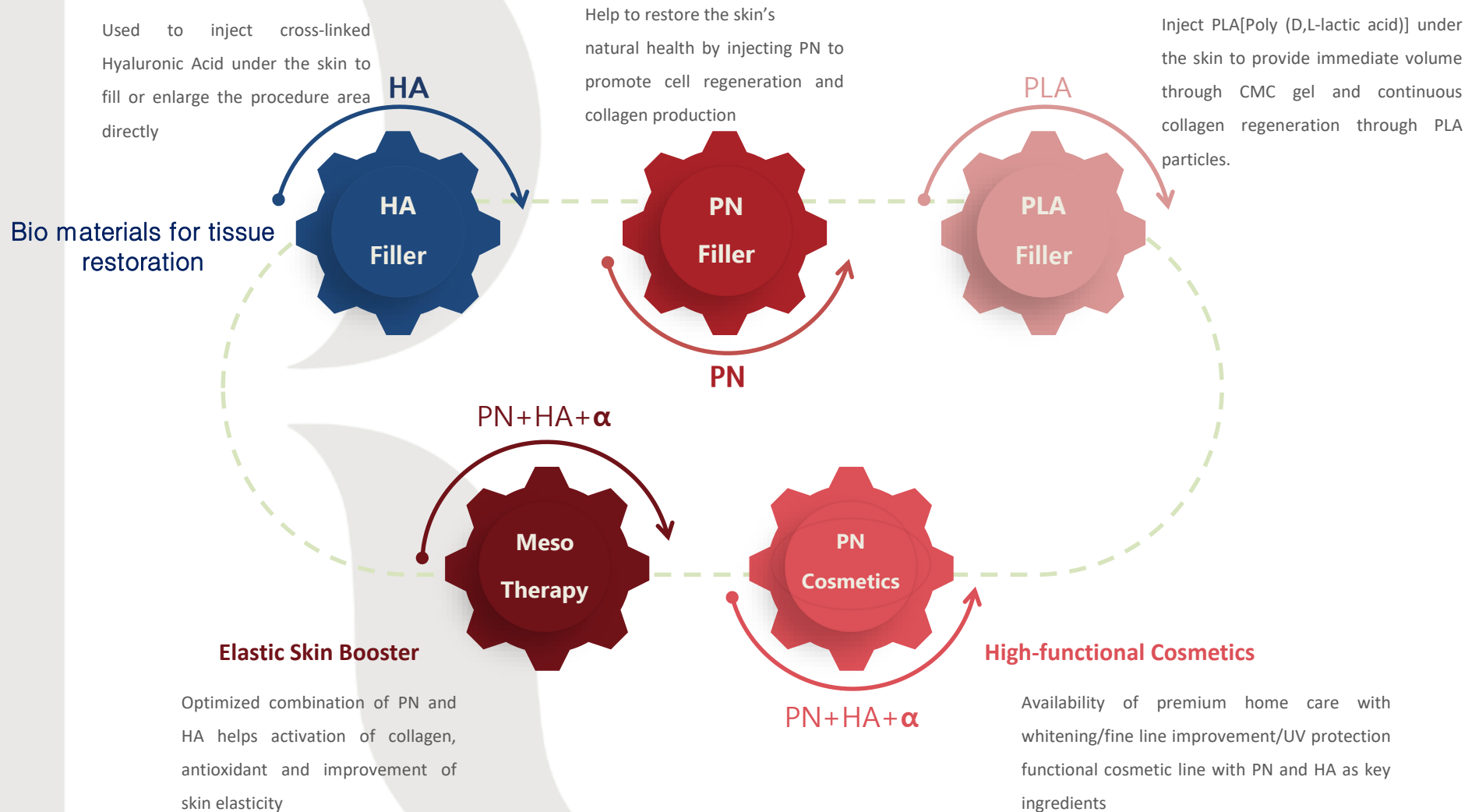

What is REVS?

A Solution line that recovers the healthy beauty of the skin by healing based on PN Ingredient extracted from salmon DNA

Special anti-aging brand with lineups covering all phases of the solution for beautiful skin, from skin booster and filler products to post-procedure soothing mask packs and home care serum.



REVS Products Line-up

REVS YOUTH LABORATORY

Clinic Care



RMT 140HPn

Elasticity Enhancing skin booster created by the best combination of HA and PN

Home Care



Booster Serum

Skin Booster combined with Oligopeptide-5, new concept anti-aging ingredient, and HA

After Care



Hydrating Gel Mask

Moisture soothing mask that covers the sensitive skin after procedure without any gap between mask and skin to help absorb effective components

- ★ More effective when mask and serum used together!
- ★ Add cream in case of dryness



Repairing Serum

Whitening and fine lines improvement functional home care serum containing PN ingredient which helps the skin regeneration



Sparkling Mask

Home Spa Mask where carbonated bubbles wakeup pores and deliver nutrients to the skin

Improve the Quality of Human life

REANZEN Co., Ltd. has a mission to become a pioneer that contributes to the healthy beauty of the world.

We provide innovative technology beauty solutions based on biotechnology.

We are committed to contributing to the happiness and healthy life of mankind with a focus on the medical aesthetic business and the derma care business.

Optimal combination of Oligopeptide-5 and HA

REANZEN developed the 'REVS Booster Serum' by finding the optimal combination of [hyaluronic acid] and [Oligopeptide-5] that attracts moisture with R & D know-how as a bio-convergence technology company.

Therefore, we have introduced **REVS Booster Serum** at reasonable prices that can be used more safely for customers who dream of bright and elastic skin.



Oligopeptide-5 development story

AIMP1 is a protein component in the human body found in the research process of new drug development.

Experts have found that AIMP1 plays a major role in wound healing, skin regeneration, and whitening, and research has developed Oligopeptide-5, a peptide derived from AIMP1. Oligopeptide-5 stimulates fibroblasts, the 'collagen factory' in the skin, increasing collagen production and making the skin brighter and more elastic.

Oligopeptide-5, an effective anti-aging ingredient that fills collagen in the skin

Collagen, the main substance that makes up the skin, begins to decrease from the early 20s and decreases to less than half when it reaches 40s, resulting in a noticeable decrease in skin elasticity. Since there is a limit to how collagen is applied to the skin and ingested, Oligopeptide-5 is the ingredient developed to supplement it. Oligopeptide-5 restores the function of fibroblasts that make collagen in the skin and helps activate collagen production.





✓ **Collagen production** by activating fibroblasts

✓ **Skin regeneration** derived from raw anti-aging ingredients and **homeostatic substances**

Growth of fibroblast collagen and KGF (keratinocyte growth factor)

As a result of measuring the collagen growth effect in KGF, it can be seen that the growth effect is better than that of TGFb1.

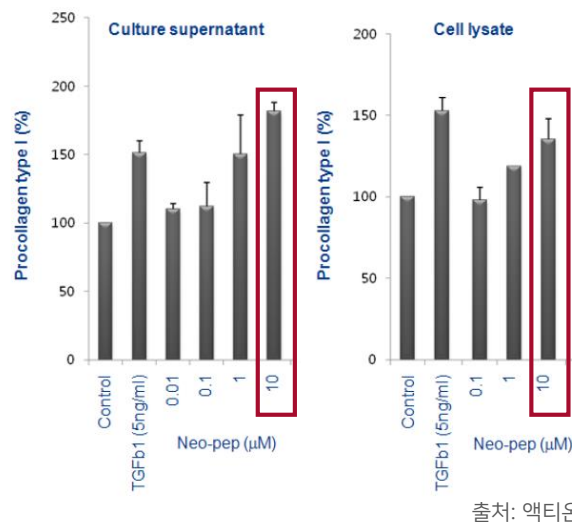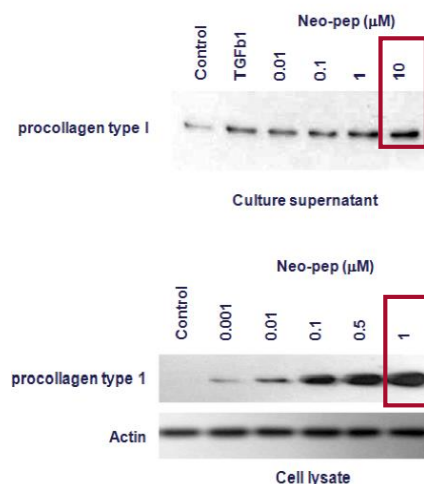


* TGFb1: Transforming growth factor beta 1,
Factors that stimulate growth factors to regenerate
and differentiate tissues

Measurement of collagen production effect compared to EGF

Oligopeptide-5 can care for skin aging (wrinkle, sagging, complexion) due to the reduction of fibroblasts, and it is possible to improve skin aging even with a small amount compared to EGF.



In vitro gene level (Hs68, Human dermal fibroblasts)
Evaluation against Human EGF(2 ppm)

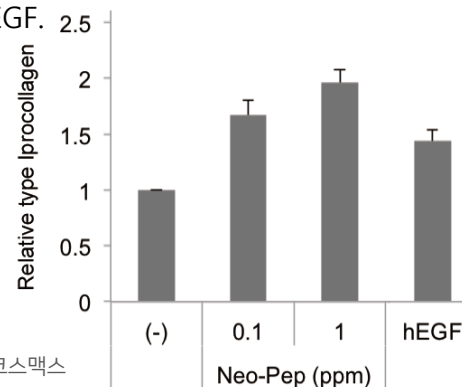

Collagen expression (candidate/hEGF)		
hEGF	Neo-Pep (ppm)	
	0.1	1
1.00	1.16	1.37

REVS' PN

The PN component in Booster Serum was developed by Reanzen's proprietary technology and is more effective than PDRN. REVS Booster Serum maximizes skin regeneration ability, restores elastic fibers by aging, and maintains skin moisture balance.

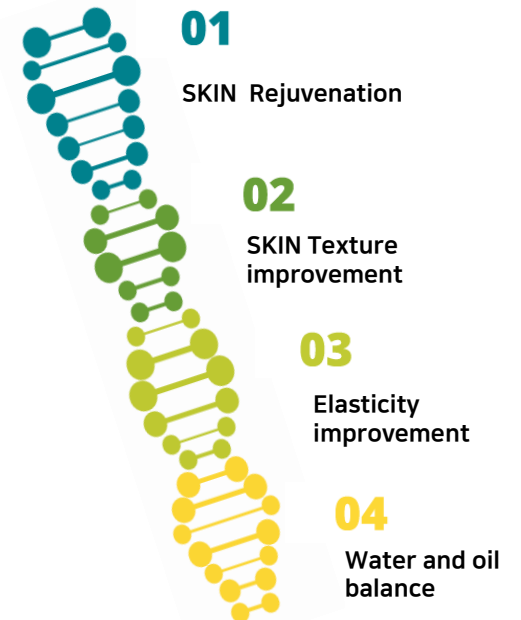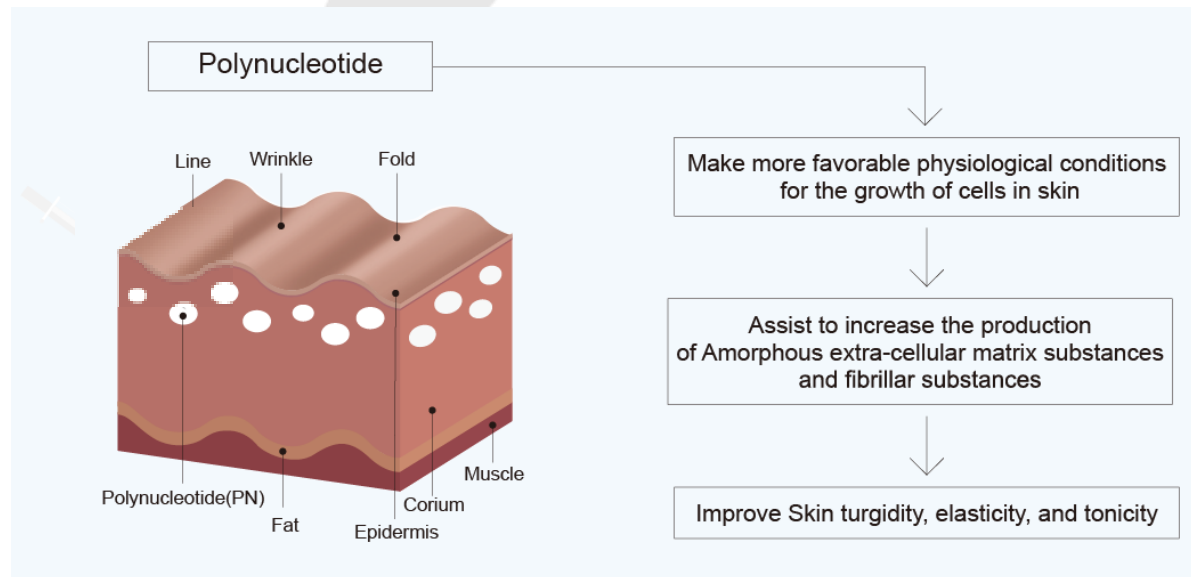

Why REVS Booster Serum?

REVS Booster Serum fills collagen in the skin.

Collagen begins to decrease from the early 20s and decreases to less than half by the 40s, resulting in a noticeable decrease in skin elasticity. The most effective resilience care method is to restore the function of fibroblasts that make collagen in the skin and help activate collagen production.

REVS Booster Serum works through the REVS Connect Roller to work on the skin epidermis, helping to fill the collagen in the skin more effectively than the commonly applied cosmetics.



REVS Booster Serum 1set(4mL*5ea)

REVS Booster Serum is a skin booster product for spa / skin care shops.

Oligopeptide-5, which stimulates fibroblasts that produce collagen in the skin and helps skin elasticity, is the main ingredient of Booster Serum, and it is an optimal combination with hyaluronic acid to help improve skin tone and elasticity. The greatest effect is achieved when treated with REVS Connect Roller.



REVS Booster Serum Device

REVS Microdisc Needle Care is a device product for mesotherapy that can effectively improve skin absorption by using a disc roller with about 210 fine needles to aid skin penetration.



REVS Connect Roller

REVS Connect Roller




Easy Use

Vial connecting design that can be easily used during home care without inconvenience of applying ampoules separately



Safe Disc Needle

Safe needle specialized for home care in the form of a disc needle in which the pain-minimized triangular needle rotates individually

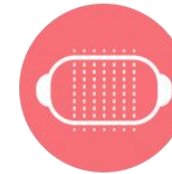

Micro-Needle Care

Healthy Skin Irritation with a length of 0.2mm and 210 micro-needles



Gamma Sterilization completed

Double sterilized and safely packaged with clean use

※The yellow part is a color change due to gamma-ray sterilization, so you can use it with no worries

REVS MTS Connect Roller delivers REVS' specialized ingredients to the skin

REVS Connect Roller

With patented technology **LOW** pain **HIGH** effect!



Triangular Needle
Smaller Fine holes,
Minimize pain and
dead skin cell



Disc Needle
Minimize pain and
wounds with individual
spinning disks




One-touch Fluid-control device
Easily adjust amount of
solution with light touch



Flow Path Disc Needle
As the disc rotates, the
ampoule is delivered by
centrifugal force through
the flow path of the
needle surface

Product with 30 individual patents

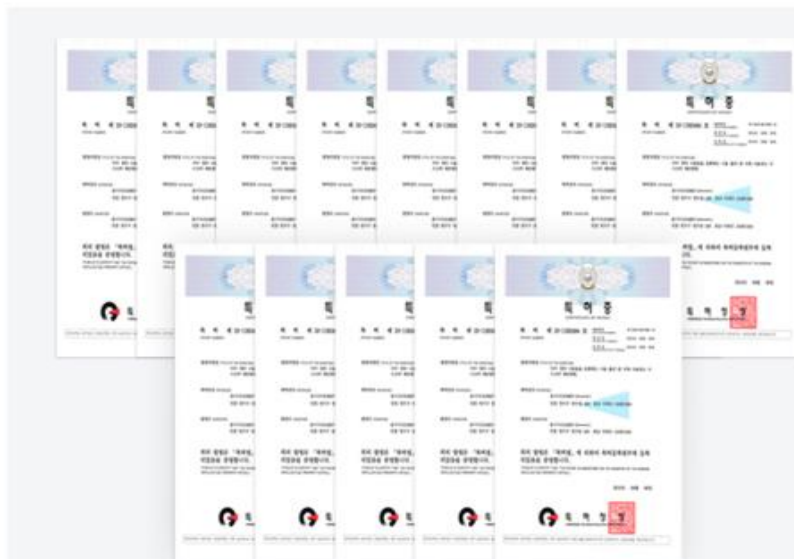

Patent on manufacturing method of triangular needle(No. 10-1182494)

Skin Stimulator Patent(No. 10-1383203)

Roller Assembly Patent(No. 10-1686226)

Beauty equipment related Patent(No. 10-1647009)

Micro needle Patent(No. 10-1112095)

Micro Needle Roller Patent(No. 10-09088893)

Multiple Needle Roller Patent(No. 10-1120467)

Patent for Fluid Applier(No. 10-1646934)

Beauty Equipment Patent(No. 10-1647836)



Patent for connection caps for needle equipment(No. 10-1658229)

Patent for needle roller using fluid(No. 10-2016108)


Patent for Cushion Needle Roller Head(No. 10-1108992)

Patent for Micro Needle Roller (No. 10-0959075)


REVS Connect Roller Principles of Action

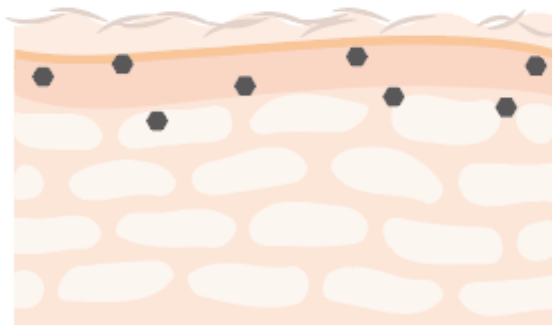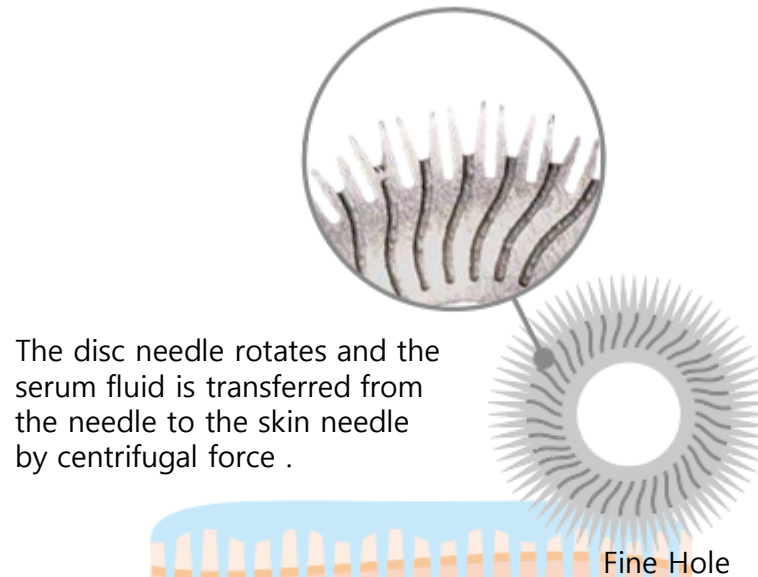

Skin trouble / wrinkles / rough skin

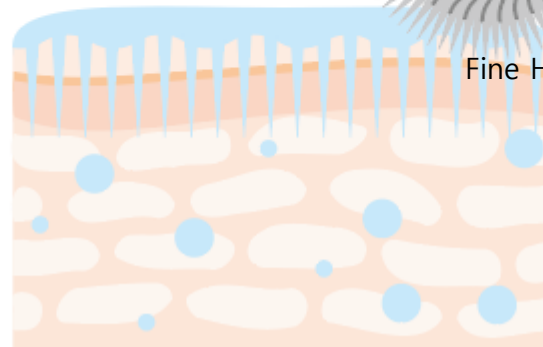

Absorb active ingredients

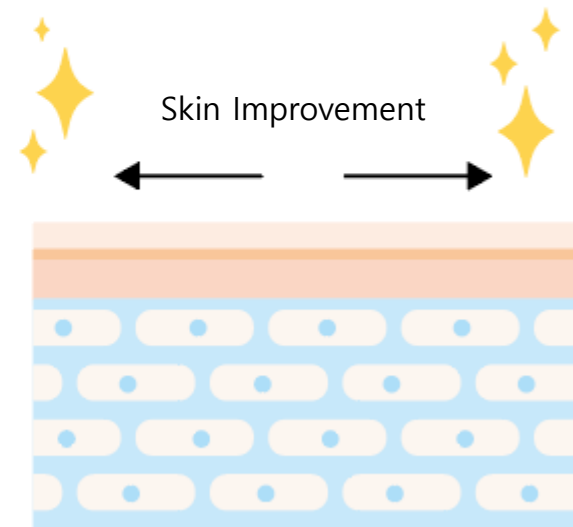

Promote Collagen Production

Skin Before using REVS >

REVS Booster Serum >
+ Connect Roller

Skin after using REVS

How to use REVS Connect Roller



1. Safely remove the lid of the REVS Booster Serum



2. Lightly put on the bottle into the silicon insert of the REVS Connect Roller



3. Make sure that the joint is connected without any problems and open the lid to prepare it



4. Tap the bottom of the container lightly with your index finger and adjust the amount of solution.

REVS Booster Serum Usage Cycle

Use **REVS Booster Serum**
once a week to experience a
noticeable change in skin
over a total of five weeks

**EVREY 7DAYS,
REVS Booster Serum
1 Vial**



Thank You

[REANZEN CO., LTD]

T. 031-596-6310

F. 031-596-6311

www.reanzen.com