



MDP CO., LTD.

19F, Gasan Digital 1-ro 88, Geumcheon-gu, Seoul, KOREA (08590)

Tel +82-2-6419-7007 Fax +82-2-6419-7015

<http://www.bravomdp.com>



Smart, Secret, Simple for Body care

LaBD LU





We change the myth that diet is difficult just for you. Lovely women, LU, a Beauty health care device, will be with you wherever whenever. LU makes you Lovely You, always wanting to be beautiful and loved. LU, the Beauty master, helps you manage your abdomen and other body parts.




LaBD LU

LaBD LU is a Home-Use Aesthetics Device using Low Level Laser helping customers manage their abdomen areas and overall body shape. Unlike traditional analogue laser type, LaBD LU lasts semi-permanently and does not need any consumables and its volume is minimized. Low Level Laser delivers the stable and exact laser power and can be worn in certain body parts without damaging the skin texture.



LaBD LU How to Use

When to use

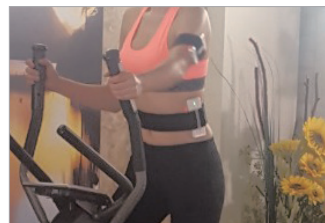

① Before and After work-out (Good for pre-exercise warm up and enhance exercise effect) ② Wearable in the work place, watching TV or Movie ③ Wearable while cleaning or dish washing at home ④ On your body, anywhere you want to manage. (Abdomen, Flank, Arms, Thighs, Chin)

On your body, anywhere you want to manage. (Abdomen, Flank, Arms, Thighs, Chin)

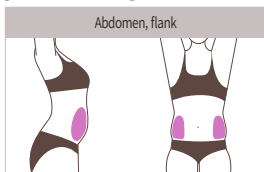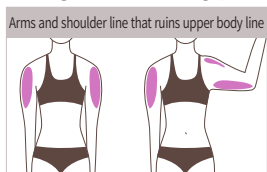


- Helps reduce skin surface roughness
- Helps tighten problem area(sagging arms, thighs, stomach)
- When used during work out, it helps facilitate metabolism and reduce edema




LaBD LU EMC(Electromagnetic Compatibility)




KC

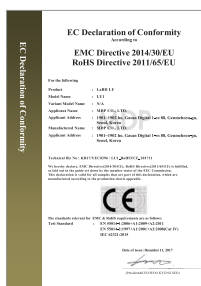


CE

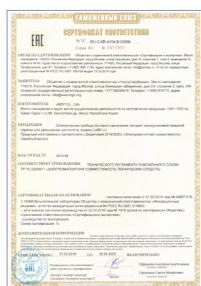


CU

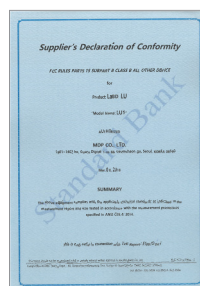

FCC

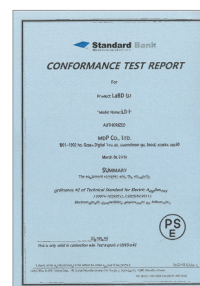

PSE

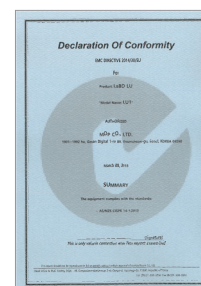

C-Tick

World widely EMC certified (Korea KC, EU CE, US FCC, Japan PSE, Australia C-tick and Russia 3 nations CU) and safe product!

LaBD LU Product Specifications & Performance

- Laser Class : Class 3R (Low-level laser)
- Laser Diode Power output : 5mW(Max)
- PulseWidth/Frequency : 50 us/4.46kHz
- Laser wavelength : 650~658nm
- Battery specification: Lithium-ion battery, 3.7V / 1700mAH
- Dimension: 120mm*50mm*130mm
- Weight: 130g
- Characteristics: Portable & Re-chargeable
- Materials : ABS
- Colors : Silver

- Operation temperature : 0~40 C
- Charging time : About 3~4 hrs + (Depending on battery use environment)
- Charging Cable : 5V 2A Micro 5pin
- Operation mode
 - 1)Continuous operation mode
 - 2)Direct contact operation mode
- Safety device: Contact sensor
- Irradiation time: Max 30minutes (set in units of 10 minutes)

- Manufacturer : MDP Co.,Ltd
- Product name : LaBD LU
- Read the directions for "Notes" and "How to Use" carefully.