



HOMEMEDICARE SPOT PATCH  
HOMEMEDICARE SOS CLEAR SPOT PAD


homemedicare<sup>+</sup>

# BRAND STORY

Homemedicare, the beginning of smart trouble-care!

Homemedicare's perfect control from the stages from raw material preparation to packaging

Optimized solution for trouble skin Homemedicare develops the best product and proposes to you, knowing anyone's concern

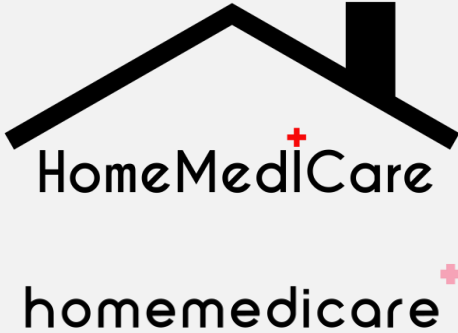

**HOME = Family(家)**  
**MEDI = Beauty(美)**  
**CARE = Health(康)**


**- A brand responsible for skin health and beauty at home-**

Homemedicare was founded in May 2016 and began selling acne patches targeting the Chinese market.

So far, we have been steadily expanding into overseas markets, including China, Taiwan, Hong Kong, Thailand, the United States and the UK. Especially, the Homemedicare Spot patch products have been loved by consumers. We had a great response as our product recorded complete sold-out in kakao Makers in Korea. In addition, we won the idea-product (a novelty) certification in the 2018 Hi-Seoul Excellent Product Brand Award, thus getting the quality of our Homemedicare products recognized.


In the future, Homemedicare will keep launching new products with differentiated quality by conducting meticulous control from the stages from raw material preparation to packaging

A character that can imprint the identity  
 (friendly and loving images) of our brand on consumers  
 <Mandu character - KANGKANG / MEIMEI>  
 Strengthening our brand image through the character

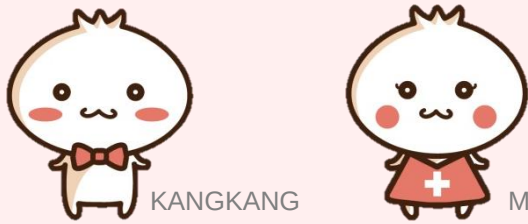

KANGKANG      MEIMEI

## CERTIFICATE & PATENT & AWARD
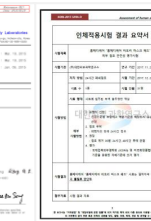
|                                                                                                                              |                                                                                                                                                   |                                                                                                             |                                                                                                             |                                                                                                                                         |                                                                                                                                                    |                                                                                                                                                 |                                                                                                                       |
|------------------------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------|
|  <p>HERBFCF Patent<br/>No. 10-1571323</p> |  <p>HomeMediCare Spot<br/>Patch CFDA<br/>Sanitation Permit</p> |  <p>FDA registration</p> |  <p>KGMP certificate</p> |  <p>Homemedicare<br/>Spot Patch<br/>Certificate</p> |  <p>Homemedicare<br/>SOS Clear Spot Patch<br/>Certificate</p> |  <p>CGMP</p>                                               |                                                                                                                       |
|  <p>ISO 22716</p>                         |  <p>ISO 14001</p>                                              |  <p>ISO 13485</p>        |  <p>ISO 9001</p>         |  <p>SPF Test<br/>Report</p>                          |  <p>Skin Sensitization<br/>Test Data</p>                       |  <p>Cytotoxicity,<br/>Skin Sensitization<br/>Test Data</p> |  <p>Skin<br/>Irritation Test</p> |



# PRODUCTS

## HOMEMEDICARE SPOTPATCH PRODUCTS

- Product Name: Homemedicare Spot Patch (Acne Patch)
- Capacity: 48 sheets
- Expiry date: 36 months from date of manufacture
- How to use: Remove the release coat according to a wound area and attach the adhesive surface of the product flat on the wound area.

**Q: What are the difference between Homemedicare Spot Patches and other companies'?**


**A:** We have complemented the three most important points that consumers care about when purchasing a trouble patch product.

1. It should be thin and not apparently recognizable.
2. It should not fall off even after you wash or shower.
3. It should be easy to carry and use and should be hygienic.



### Product Features

- Thin and not recognizable - The attached sheet must be recognizable and thin; Even after your makeup, your makeup should not be stuck into space between sheet and skin.
- Double protection and double waterproof - Hydrocolloid + pu film structure prevent the edges from getting tattered even after your wash or shower.
- Small package (domestic only) - hygienic and portable.
- Easily removable cutting film - easy to use folding.
- UV protection - SPF46 certified.
- Use of characters - Cute and lovely characters for targeting young people in their 10s and 20s.
- Reliable quality - US FDA and ISO13485 certified / passed a skin irritation test.
- \* The idea product (a novelty) certification in 2018 Hi-Seoul Excellent Product Brand Award

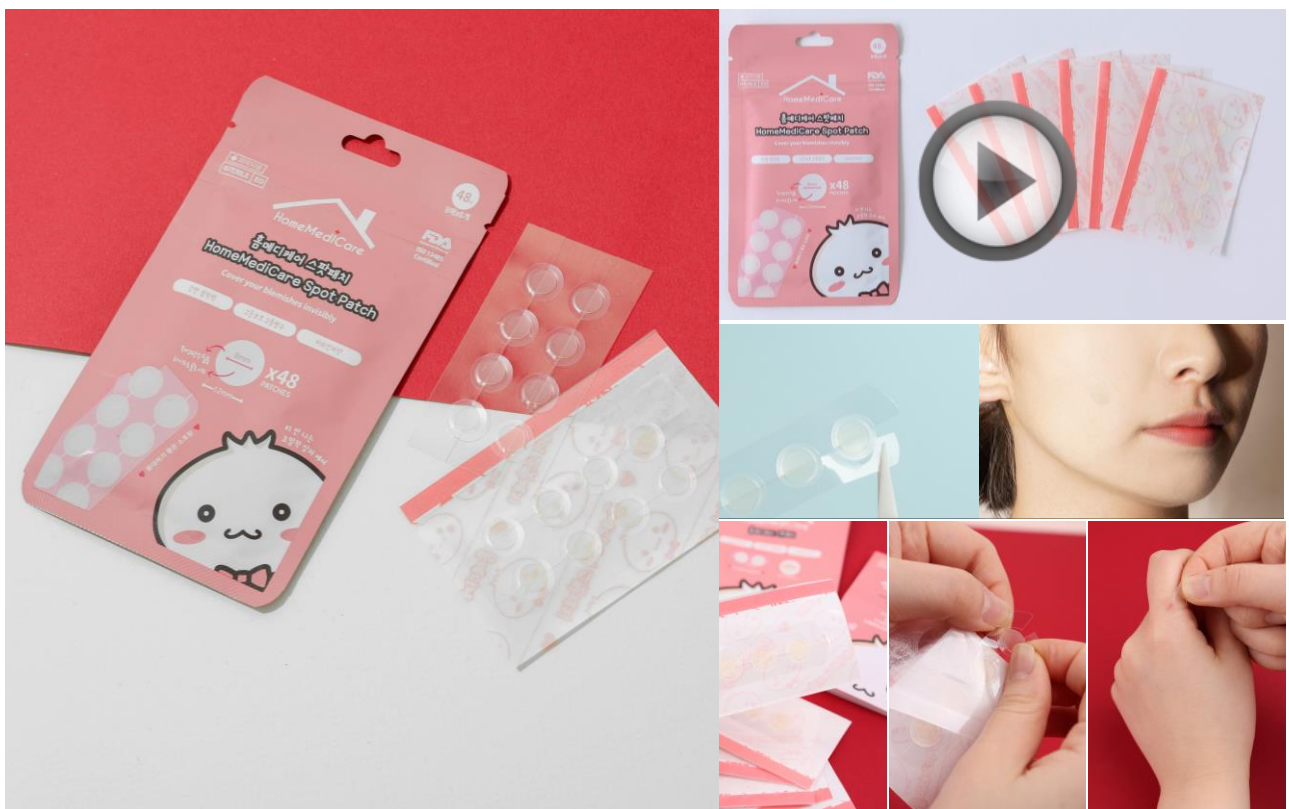

# PRODUCTS

## HOMEMEDICARE SPOTPATCH PRODUCTS

01

The edges of the attached sheet are not easily getting tattered even after your wash or shower.

The Hydrocolloid + pu film structure has a good adhesive force, preventing the sheet from falling off or getting tattered even after your wash or shower. The sheet can stay long enough to heal the wound!



Other companies'

VS



Homemedicare Spot Patch

02

Thin and not easily recognizable, making your makeup OK!

The attached sheet is not easily recognizable and thin; even after your makeup, your makeup must not be stuck into space between sheet and skin, enabling you to have a natural image



Other companies'

VS



Homemedicare Spot Patch

03

More hygienic and more convenient Homemedicare spot patch

Its individual small package enables you to always carry with you in your wallets, handbags, pouches, and phone cases! Resolve your skin concern hygienically anywhere, anytime!



"Your must-have item in a pouch! You do not need a concealer! "



Other companies'

VS



Homemedicare Spot Patch

04

Cutting film - easy to remove, no edge folding


Hygienic, easy to use! You do not have to forcibly remove it, you do not have to worry about edge folding.

" Say Goodbye to unsanitary spot Patch products~~ !! "



Other companies'

VS

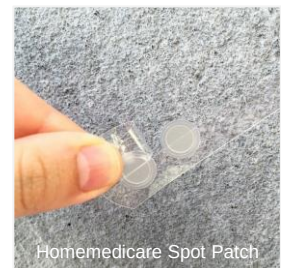

Homemedicare Spot Patch

### Completion of the test for safety and UV protection!


Giving you an SPF 46 Sunscreen effect! Solve your skin anxieties and protect your already weakened skin from ultraviolet rays.




"Because my skin is precious!"

# PRODUCTS

## HOMEMEDICARE SPOTPATCH PRODUCTS

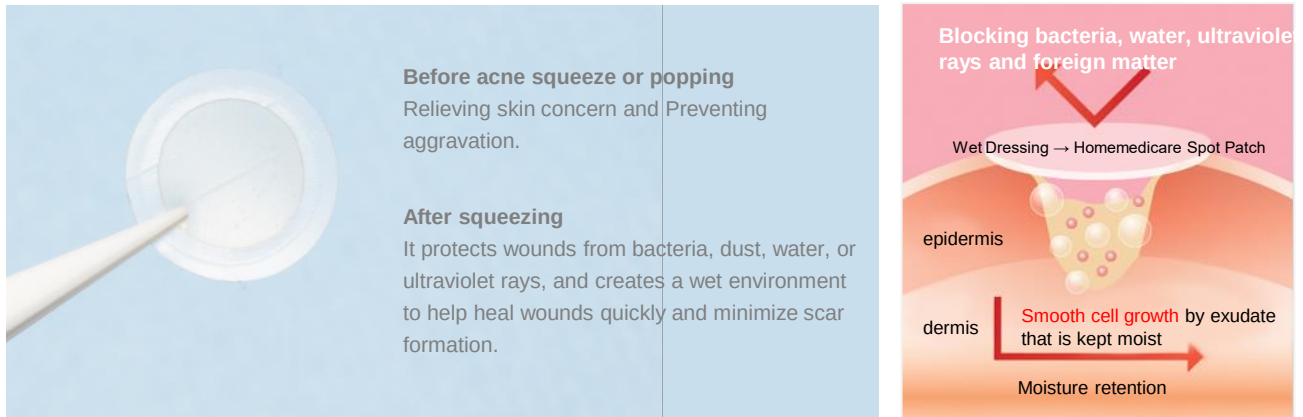

## HOME MEDICARE SPOT PATCH HEALING PROCESS

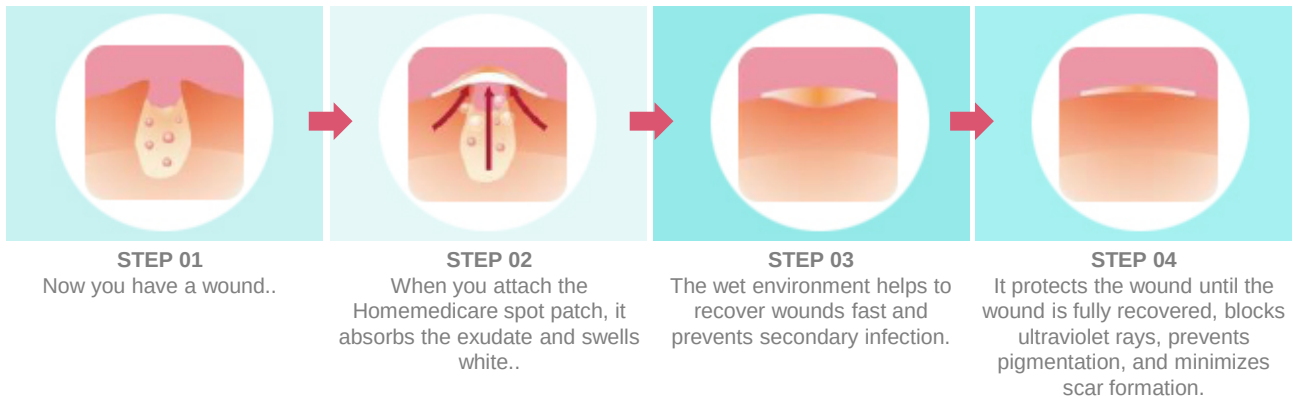

## HOW TO USE THE HOMEMEDICARE SPOT PATCH

### Replacement Period

If the product swells white, the exudate leaks out, or the whole product turns white, please replace it. (If the product is swells whitened, it means that the wet environment is well formed and the healing process is under way. Please maintain the condition as long as you can.)


### Replacement Method

When replacing, grasp the end of the product and remove it softly while trying keeping the sheet horizontal with the wound surface.

If it is removed forcibly, it may damage your regenerated skin tissue. So, please apply warm water to soften the patch and remove it.

### Precautions

1. Do not use it for infected wounds or deep wounds.
2. If you see itching, rash, or redness, discontinue the use and consult a specialist.
3. Do not reuse used products.
4. Do not apply a disinfectant (such as ointment, cream, or iodine solution) on the wound together with this product.
5. If the product needs to be removed early, remove it after applying lukewarm water and making the sheet swollen.
6. Do not use it for patients with diabetes or blood circulation troubles and those with easy-to-be-rash skin.

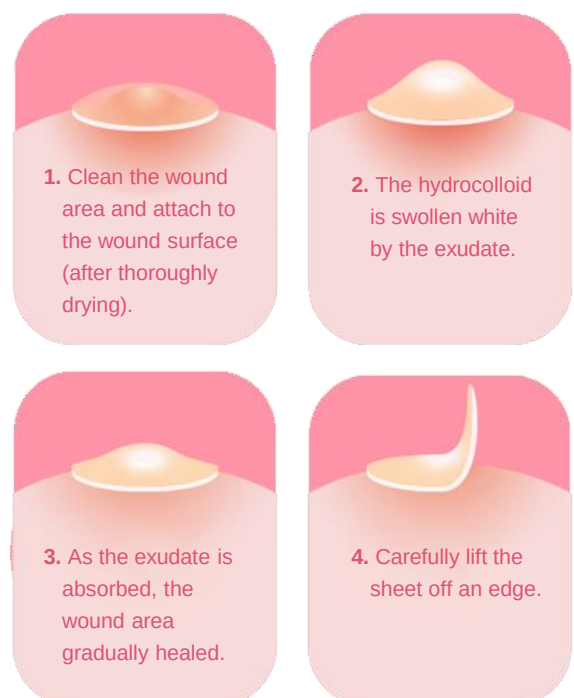

# PRODUCTS

## HOMEMEDICARE SOS CLEAR SPOT PAD

- Capacity: 150ml / 70 sheets
- Expiry date and Expiry date after opening: Expiration date is 3 years from the date of manufacture. Recommended use within 12 months after opening
- How to use: After wiping the cotton pad to a direction of skin texture, lightly tap it and make it absorbed. (Tip! When your skin is extremely sensitive and dry, it will be more effective if you put it on the sensitive area and remove after 5 ~ 10 minutes)

**Q: Why did you make a sos clear spot pad?**

**A:** Most people who are worried about skin troubles tend to have troubles on partial areas, such as cheeks, foreheads or jaw spots, not on a whole face. However, it is said that a whole-face mask sheet with strong ingredients covering a whole face is burdensome. That is why Homemedicare has created a trouble-focused intensive care product with a natural tea tree ingredient and Herb CFC (Herb CFC) which can be applied only to the part of trouble..



SOS CLEAR SPOT PAD

### Product Features

- It specializes and strengthens a natural tea tree component, which is universally the most effective ingredient for skin trouble improvement.
- Patented No. 10-1571323 <Cosmetic compositions comprising extract of herbs: synergy effect created by Proven skin care ingredient tea tree + Patented trouble care ingredient (Herb CFC) >
- All ingredients are US Skin Deep EWG Green grade material (all good ingredients)
- Using natural ingredients
- 100% cotton pad with no irritation perfectly sticks to your skin
- Completion of skin irritation test – Result of no irritation
- \* Korea Dermatology Research Institute / Number of test subjects: 31 persons / Test result: Judged as non-irritating substance.
- No artificial fragrance - natural tea tree fragrance soothes both your skin and mind
- Contains tongs for hygienic use
- Trustworthy quality competitiveness - Manufactured under certified technology and quality control.
- CGMP / ISO9001 / ISO14001 / ISO22716 Certification, the idea-product certification in the 2018 Hi-Seoul Excellent Product Brand Award





# PRODUCTS

## HOMEMEDICARE SOS CLEAR SPOT PAD

Refreshing your sensitive skin without stickiness! Free of artificial fragrance, natural tea tree fragrance calms both your skin and mind!

**Natural tea tree extract and 7 vegetable complex herbal extract from Jeju makes a sensitive and dry skin into a comfortable and moist skin** <Synergy effect of tea tree, a patented trouble care ingredient + certified trouble care ingredient!>



**Tea Tree**  
Disinfection effect that is 20 times as strong as alcohol, sterilization effect. Helps to calm and soothe skin.



**Peppermint**  
Cleanses pores and has excellent astringent. Helps to calm and moisturize skin



**Rosemary**  
Gives skin vitality and vitality, helps prevent moisture evaporation.



**Chamomile**  
Inhibits elastin degrading enzyme and regenerates skin. Gives skin vitality and energy, helps to prevent water evaporation



**Lemon grass**  
skin cleansing. Helps with pore convergence and sebum control.



**Mother chrysanthemum**  
excellent in moisturizing power. Inhibits the activity of substances harmful to the skin and removes them to help calm the problematic skin



**Lavender**  
skin disinfection. Helps to calm and relax trouble-prone skin



The patented vegetable complex extract helps to soothe the troubled skin.

### Good ingredient

The US Skin Deep ® EWG green-grade ingredients (EWG rating is in each bracket)

Purified water (1), Dipropylene glycol (1), Tea tree leaf (1), Betaine (1), 1,2-hexanediol (1), Tea tree leaf oil (1), polyglutamic acid (1), Polyglyceryl-10 laurate (1), Polyglyceryl-10 myristate (1), Sodium cocamidopropylphosphorylcholine phosphate (1), polyhydroxystearic acid (1), Isononyl Isononanoate (1), Ethylhexylisononanoate (1), Black Willow bark extract (1), Butylene glycol (1), Rosemary extract (1), Lavender flower extracts (1), Tea tree leaf extract (1), Chamomile flower extract (1), Wild Chrysanthemum flower extract (0), Lemon grass extract (1), Peppermint leaf extract (1), Xylityl glucoside (1), Anhydroxylitol (1), Xylitol (1), Glucose (1), Disodium di-t-iodide (1), Chlorphenesin (2)



Its 100% cotton pad makes it cleaner and smoother



Tongs for hygienic use

KOSHA 2017-10-04-01

Assessment of human skin compatibility

인체적용시험 결과 요약서

|       |                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                            |       |                             |  |  |
|-------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------|-----------------------------|--|--|
| 시험목적  | Homemedicare® ‘클리어스팟 패드’의 피부 자극성 평가<br>피부 자극성 안전성 평가하기                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                     |       |                             |  |  |
| 시험기관  | 한국소비자원                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                     | 시험 일자 | 2017. 11. 23 ~ 2018. 01. 05 |  |  |
| 시험 방법 | 인체 피부 자극성 시험                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                               | 시험 일자 | 2017. 11. 23 ~ 2017. 12. 05 |  |  |
| 시험 대상 | 건강한 성인                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                     | 시험 인원 | 31명                         |  |  |
| 시험 결과 | <div> <div>대부분의 시험대상자</div> <div> <div>1. 피부 자극성 평가 결과</div> <div>2. 피부 자극성 평가 결과</div> <div>3. 피부 자극성 평가 결과</div> <div>4. 피부 자극성 평가 결과</div> <div>5. 피부 자극성 평가 결과</div> <div>6. 피부 자극성 평가 결과</div> <div>7. 피부 자극성 평가 결과</div> <div>8. 피부 자극성 평가 결과</div> <div>9. 피부 자극성 평가 결과</div> <div>10. 피부 자극성 평가 결과</div> <div>11. 피부 자극성 평가 결과</div> <div>12. 피부 자극성 평가 결과</div> <div>13. 피부 자극성 평가 결과</div> <div>14. 피부 자극성 평가 결과</div> <div>15. 피부 자극성 평가 결과</div> <div>16. 피부 자극성 평가 결과</div> <div>17. 피부 자극성 평가 결과</div> <div>18. 피부 자극성 평가 결과</div> <div>19. 피부 자극성 평가 결과</div> <div>20. 피부 자극성 평가 결과</div> <div>21. 피부 자극성 평가 결과</div> <div>22. 피부 자극성 평가 결과</div> <div>23. 피부 자극성 평가 결과</div> <div>24. 피부 자극성 평가 결과</div> <div>25. 피부 자극성 평가 결과</div> <div>26. 피부 자극성 평가 결과</div> <div>27. 피부 자극성 평가 결과</div> <div>28. 피부 자극성 평가 결과</div> <div>29. 피부 자극성 평가 결과</div> <div>30. 피부 자극성 평가 결과</div> <div>31. 피부 자극성 평가 결과</div> </div> </div> |       |                             |  |  |
| 시험 결과 | <div> <div>대부분의 시험대상자</div> <div> <div>1. 피부 자극성 평가 결과</div> <div>2. 피부 자극성 평가 결과</div> <div>3. 피부 자극성 평가 결과</div> <div>4. 피부 자극성 평가 결과</div> <div>5. 피부 자극성 평가 결과</div> <div>6. 피부 자극성 평가 결과</div> <div>7. 피부 자극성 평가 결과</div> <div>8. 피부 자극성 평가 결과</div> <div>9. 피부 자극성 평가 결과</div> <div>10. 피부 자극성 평가 결과</div> <div>11. 피부 자극성 평가 결과</div> <div>12. 피부 자극성 평가 결과</div> <div>13. 피부 자극성 평가 결과</div> <div>14. 피부 자극성 평가 결과</div> <div>15. 피부 자극성 평가 결과</div> <div>16. 피부 자극성 평가 결과</div> <div>17. 피부 자극성 평가 결과</div> <div>18. 피부 자극성 평가 결과</div> <div>19. 피부 자극성 평가 결과</div> <div>20. 피부 자극성 평가 결과</div> <div>21. 피부 자극성 평가 결과</div> <div>22. 피부 자극성 평가 결과</div> <div>23. 피부 자극성 평가 결과</div> <div>24. 피부 자극성 평가 결과</div> <div>25. 피부 자극성 평가 결과</div> <div>26. 피부 자극성 평가 결과</div> <div>27. 피부 자극성 평가 결과</div> <div>28. 피부 자극성 평가 결과</div> <div>29. 피부 자극성 평가 결과</div> <div>30. 피부 자극성 평가 결과</div> <div>31. 피부 자극성 평가 결과</div> </div> </div> |       |                             |  |  |

Homemedicare

3

본 시험결과서와 “인체적용시험 결과보고서”는 시험대상자의 동의하에 작성된 것으로, 시험대상자의 동의 없이 시험결과서와 “인체적용시험 결과보고서”를 제3자에게 제공하지 않습니다.

### Sebum control / moisturization / calming effect without stimulation (no stimulation)

- Test period: December 26, 2016 ~ 12.29
- Number of test subjects : 31 persons
- Institution: Korea Dermatology Research Institute
- Test result: judged to be non-irritant

### Test Result

The samples of Homemedicare SOS Clear Spot Pad are judged to be non-irritating substances in terms of the primary stimulation standard.

## HOW TO USE

### STEP 01

After cleansing, it is used for the skin, toner stage. Wipe it once more with the pure cotton pad to wipe away any remaining makeup or impurities.



### STEP 02

When the skin is extremely sensitive and dry, it is more effective to put it on the sensitive area and remove it after 5 ~ 10 minutes to let it sufficiently absorb the product





**HomeMediCare Co.**

A-603, 45, Magokjungang-6ro, Gangseo-gu, Seoul, Republic of Korea

T. +82-70-8624-9033      E. homemedicare@naver.com

homemedicare.co.kr