

Contents

01 Product Features_Handpiece

02 Product Features_Disposable Tip

03 Features of the Company

01

Product Features_Handpiece



: Advantage in 'Technology' and 'Design'



Technology



Design








- 1) Technology_Powerful 'BLDC Motor(Brushless DC Motor)'
 - Low mechanical vibration, low noise even at high speed
 - Durability that can be used for more than 3 years
 - Having more than 4 times stronger power than general motor, it can be used without skin damage in case of deep skin treatment

	AC motor	DC motor	BLDC motor
Durability	High	Low	Semi-permanent
Noise	Noisy	Noisy	Low
Direction of rotation	Bidirectional	Unidirectional	Bidirectional
Power consumption	High	Low	Low

2) Technology_Stable Power Supply

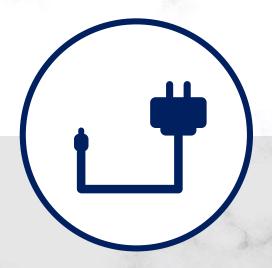

<Wired type>

Stable power delivery during treatment



<100V~250V Free volatage>

Can be used in any country



<Length of cable:4.4m>

Can be used freely
with long cable adapter
(extention cable: 1.8mm / adaptor: 2.6m)




3) Design_Ergonomic Design



- Precise needle length adjustment: the needle length can be adjusted in 0.1mm increments.
- Anti-slip design : Rubber cover prevents slipping during treatment
- Easy storage: Handpiece holder is provided for easy and hygienic storage after use
- Golden ratio in design: Moderate weight and comfortable grip with center of gravity at the bottom of the handpiece
- Dial type of speed adjustment: Easily adjustable while in use



- 4) Strict Quality Control
 - Product test is conducted based on actual use in a very harsh environment

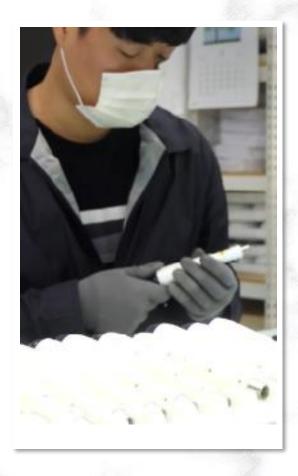









02 Product Features_Disposable Tips(Fastening Method)

1) 'Double Locking System' Using Double Threaded Screw




- Double locking fastening system
 - : Prevents needle loosening and needle popping up by vibration
- Immediately delivers motor force to the needle and adjusts needle movement from the motor force
 - : Using robust fastening method, enables skin not to be damaged or make a scratch

02 Product Features_Disposable Tips(Design)

2) Design_Patented Crown Type





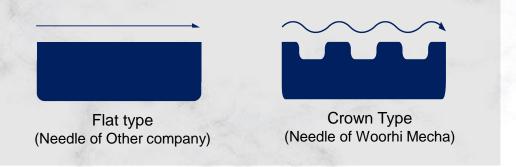


- As a special design, the area where the face and needle cartridge touch is designed as a crown type rather than a flat type.
- Prevents cosmetics / medicines from being pushed out of face during use. Free and smooth movement when needling, effective and easy to use



02 Product Features_Disposable Tips(Design)

3) Design_Spring Type and One-Way Design

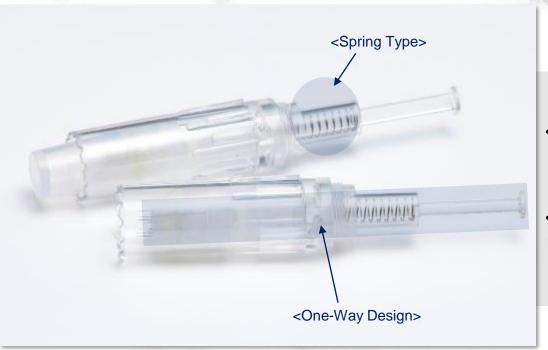

One-way design

Design that the entire spring and needle acts as an air-hole without using silicone

With air-hole, it makes air-circulation and smooth needle movement and no backflow
 The power of the motor is transferred immediately, allowing for no damage skin



02 Product Features_Disposable Tips(Design)

4) Product Assembly and Thorough Quality Control

Skilled assembler and Material management by insourcing of raw materials

- Made all processes in-house in Korea (Made in Korea)
- Own quality standard length tolerance of 0.05mm
- No needle hook / bending
- · No foreign material in the product
- Relatively price advantage in high quality production standard











02 Product Features_Disposable Tips(Quality)

4) Product Assembly and Thorough Quality Control





• ISO13485, GMP System Management





Strict Quality Control

: Follow the regulation of AQL(Acceptable Quality Limit, S-4/Special Inspection level)



02 Product Features_Disposable Tips(Quality)

5) Stability of Raw Materials



· Needle: SS304

: It is widely used in medical instruments, has strong corrosion resistance, and can prevent metal allergy to a moderate level.



ROHS certified raw materials

: Use of safe raw materials that do not contain harmful substances (lead, cadmium, mercury, hexavalent chromium, phthalate, bromine-based flame retardants)



02 Product Features_Disposable Tips(Quality)

5) Stability of Raw Materials











- Raw material supply
 - : Produce and procure all materials in Korea

· EO gas sterilization

: EO gas with excellent sterilization ability, can be sterilized and stored for a long time without product quality change or fading



02 Product Features_Disposable Tips

6) User Friendly Design



Minimization of pain
 Sharpening angle of the needle pin reduces the pain due to the 'R' type




Anti-slip design

Non-slip design in the needle guide to prevent slipping when engaging or releasing the needle with the handpiece



02 Product Features_Disposable Tips

7) Backflow prevention (optional)

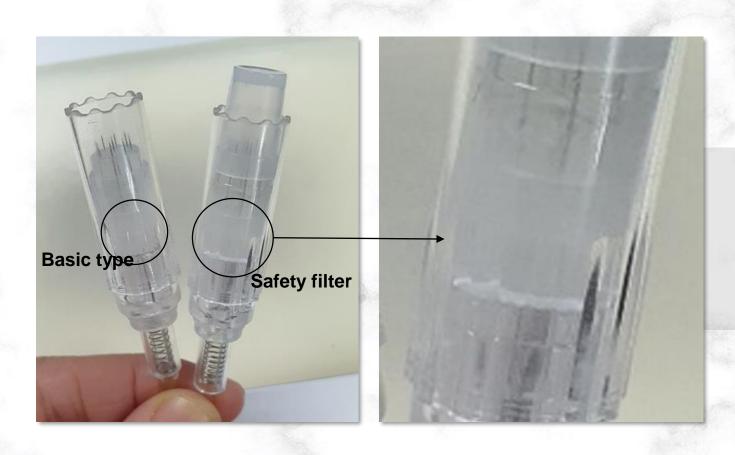

 One-way design with the entire needle does not flow back unless the hand-piece is turned upside down, but optional 'Safety filter' type for specific users



03

Features of Company



03 Product Feature_Company



- 1) Various Certificates for Company and Products
 - Certification: GMP, ISO13485, ISO 9001, ISO14001, KFDA, CE, The Best Value Certified by Gyeonggi Province
 - Awards: Prime Minister Award for Small and Medium Business Technology Innovation Competition, Award of Export Tower of 1 Million, Certificate of Merit to Creating Employment



03 Product Feature_Company

2) Overseas Distributors and Branch in Shanghai





· Overseas distributors

: We have exclusive distributors in 13 countries and have exported our products to 72 countries



Branch in Shanghai

: Provide a variety of service to customers such as distribution, exhibition, training, and after-sales service

03 Product Feature_Company



3) Having Know-how for The Product and Professional Manpower on Export

Provide know-how for product sales and procedures:
: As a manufacturer, we provide guide for product sales and effective procedures for each products





· Professional manpower on export

: we have employees who have cerificate for export business and can speak a variety of languages(English, Chinese and Spanish) for smooth communication

Thank You

